



SYNTHECON

3D Culture Systems for a 3D World

RCCS-8DQ



8 Station Rotator Base

- Rotates 8 vessels at 4 independent speeds.


Power Supply

(with Tachometer and Ribbon Cable)

Operation Manual

8 Disposable Culture Vessels

- Low cost, easy to use system which makes the technology available to most laboratories.
- Each vessel is independently packaged and gamma sterilized. Vessels are supplied 4 to a pack and are available in 10mL and 50mL sizes.

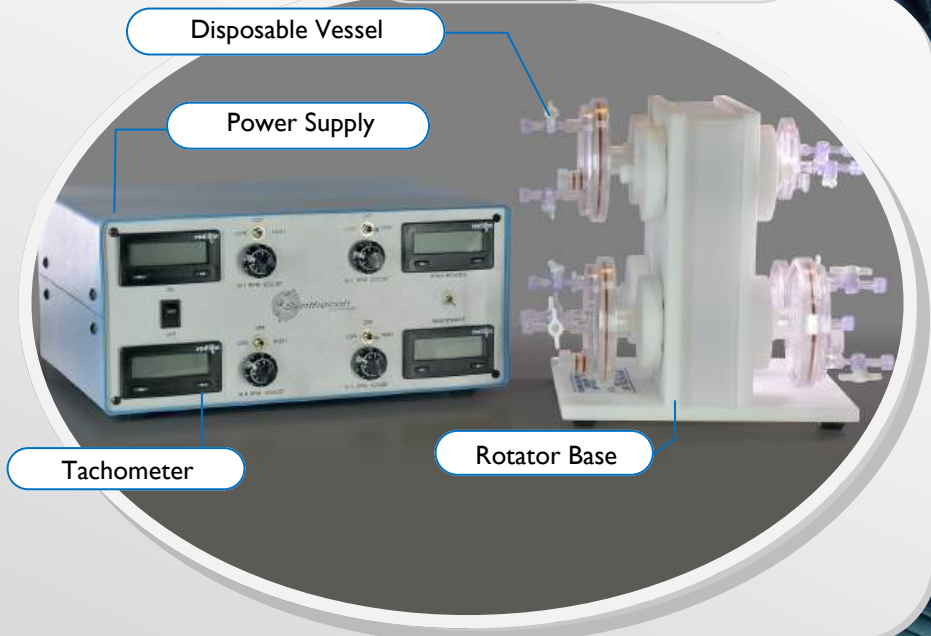

FAQ

Q: How are the cells oxygenated?

A: In a gas regulated incubator, gas diffuses through the silicone membrane of the bioreactor.

Q: Do the cells stay in one place in the bioreactor?

A: No, they continuously fall through the media during vessel rotation. The continuous motion of the cells in media facilitates their exposure to nutrients.

Q: Is the RCCS a roller bottle system?

A: No. Roller bottles grow cells on the wall of the bottle in 2D. The RCCS grows cells in suspension, either as 3D spheroids or on scaffolds/beads.

The Synthecon Rotary Cell Culture System

Produces differentiated high density, three dimensional tissue cultures.

Membrane oxygenation provides bubble free operation, which results in an extremely gentle culture.

Operates in standard incubators. Samples may be easily drawn at any time.

Ideal for growing and maintaining normal tissues, cancer tumors spheroids and bacteria/virally infected tissue models.

Synthecon can design and build systems to user's specifications.

* For additional FAQs visit www.synthecon.com/faq

RCCS-8DQ

RCCS-8DQ

TECHNICAL SPECIFICATIONS-DIMENSIONS RCCS-8DQ

COMPONENT	HEIGHT	WIDTH	DEPTH
Base	9.0" (22.9 cm)	8.0" (20.3 cm)	8.25" (21.0 cm)
Power Supply	6.0" (15.25 cm)	12.5" (31.75 cm)	9.5" (24.1 cm)
COMPONENT	WEIGHT	ROTATIONAL SPEED	—
Base	8.0 lbs (3.77 Kg)	5-35 RPM	—
Power Supply	8.0 lbs (3.77 Kg)	—	—


RCCS-8DQ POWER AND ENVIRONMENTAL SPECIFICATIONS

POWER	VOLTAGE	FREQUENCY
Power Supply Input	85-260 VAC @ .5Amps	50-60 Hz
Power Supply Output	7.5 VDC @ 2.0Amps	50-60 Hz
ENVIRONMENT	TEMPERATURE	HUMIDITY
Base	50-105 °F (10-40 °C)	0-100% RH
Power Supply	50-80 °F (10-27 °C)	0-75% RH

RCCS-8DQ COMPATIBLE CULTURE VESSELS

CULTURE VESSEL MODEL	SIZE	COMPATIBILITY
Disposable HARV *	10 mL	Yes
Disposable HARV *	50mL	Yes
Autoclavable HARV	1 mL	No
Autoclavable HARV	2 mL	No
Autoclavable HARV	4 mL	No
Autoclavable HARV	10 mL	No
Autoclavable HARV	50 mL	No
Autoclavable STLV	55 mL	No
Autoclavable STLV	110 mL	No
Autoclavable STLV	250 mL	No
Autoclavable STLV	500 mL	No

* Disposable vessels are available in the High Aspect Ratio Vessels (HARVs). Each vessel is individually packaged and gamma sterilized. Vessels are supplied four to a pack. Eight disposable HARVs are provided with the RCCS-8DQ System. All other compatible vessel types are ordered separately. Visit www.synthecon.com for more details.

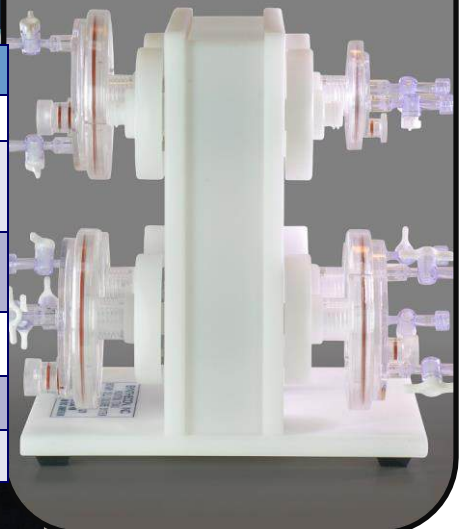

3D Culture Systems
for a 3D World

CellD SARL
19 avenue de l'Aspre
30150 ROQUEMAURE
France

0033 4 66 82 82 60

www.celld.com

contact-info@celld.com

